




MOBILE
mim™

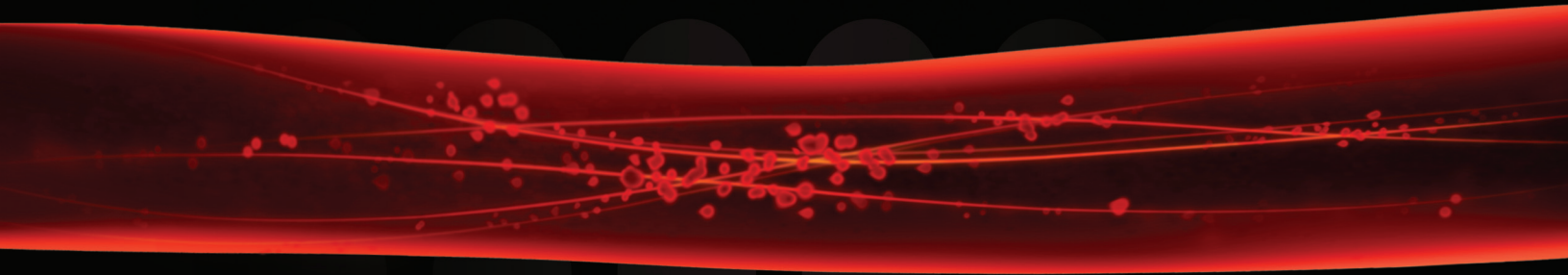

MOBILE mim™ Imaging Anywhere

Open up new possibilities with Mobile MIM.

- Mobile MIM is a diagnostic imaging app for the iPad®, iPhone®, and iPod touch®.
- Imaging professionals now have a portable solution away from their workstations.
- Consult with peers for difficult cases and overreads.
- Reduce image distribution delays.
- Enhance referring physician and patient interaction.
- Never rely on ad hoc mobile viewing solutions again.
- Bypass the cost and headache of CDs.
- Give hands-on image access to tumor board or class members.




Access your PACS with a MIM server solution.

- A MIMpacs™ server can function as an imaging server for Mobile MIM.
- Configure MIM software to cache a subset of PACS data for mobile access.
- Access patient studies with fast searching and viewing.
- Keep data within your local network.
- Access data off-site with VPN or other network extensions.

Simplify Using MIMpacs

With MIMpacs as your primary image storage solution, the setup is even simpler. Every study in MIMpacs is immediately accessible.

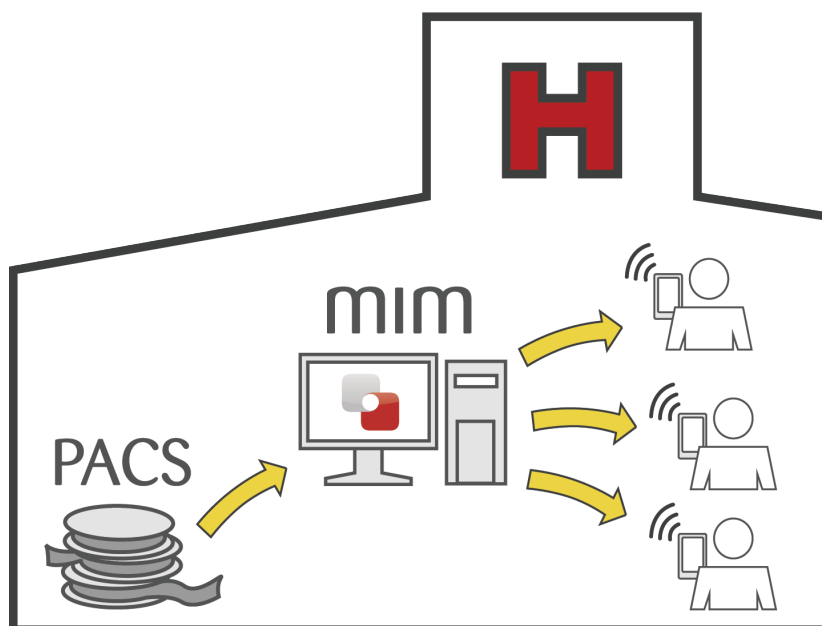


Imaging Everywhere

Extend your PACS with MIMcloud™ and MIMcloud Express™.

- MIMcloud securely hosts image data so users can maintain access while off-site.
- The free MIMcloud Express automatically routes PACS data transfers to MIMcloud.
- MIMcloud Express accepts standard DICOM transfers for institution-wide compatibility.
- MIMcloud Express is vendor-neutral and compatible with PC and Mac.
- A single server can route to multiple MIMcloud users or groups.
- All PHI is encrypted for transfer and storage (for HIPAA compliance).


MIM 5 Integrates with MIMcloud

MIMcloud compatibility is built-in to every MIM 5 workstation, allowing you to manage studies, push and pull images, and easily share studies with others inside or outside your institution.

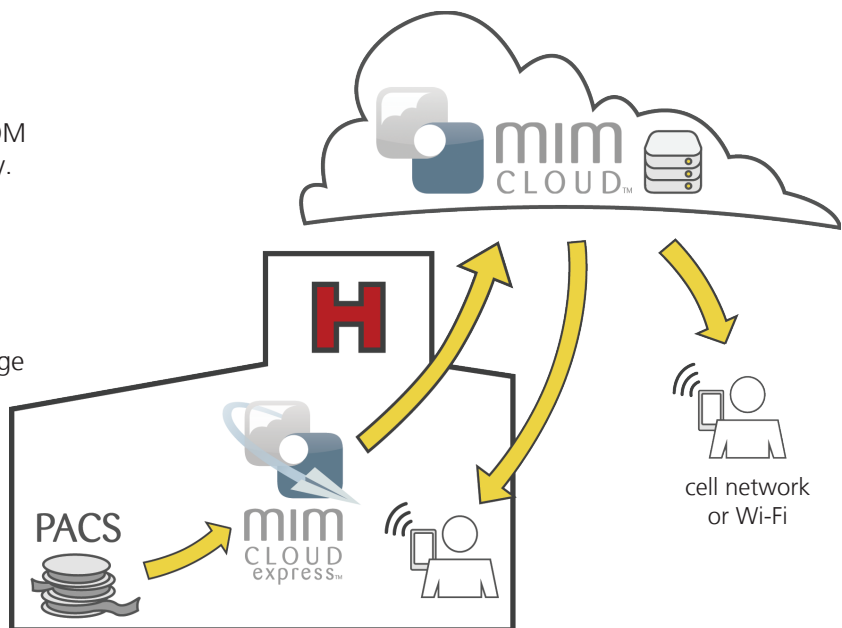

Broaden your horizons with MIMcloud.

- Access images from any Internet connection.
- Enhance teleradiology and multi-institution reading operations.
- Use MIMcloud as off-site backup.
- Share images with referring physicians, partner institutions, and patients.
- Collect or contribute images for clinical trials.
- Distribute images to users who don't have Mobile MIM using the MIMviewer® downloadable desktop software.

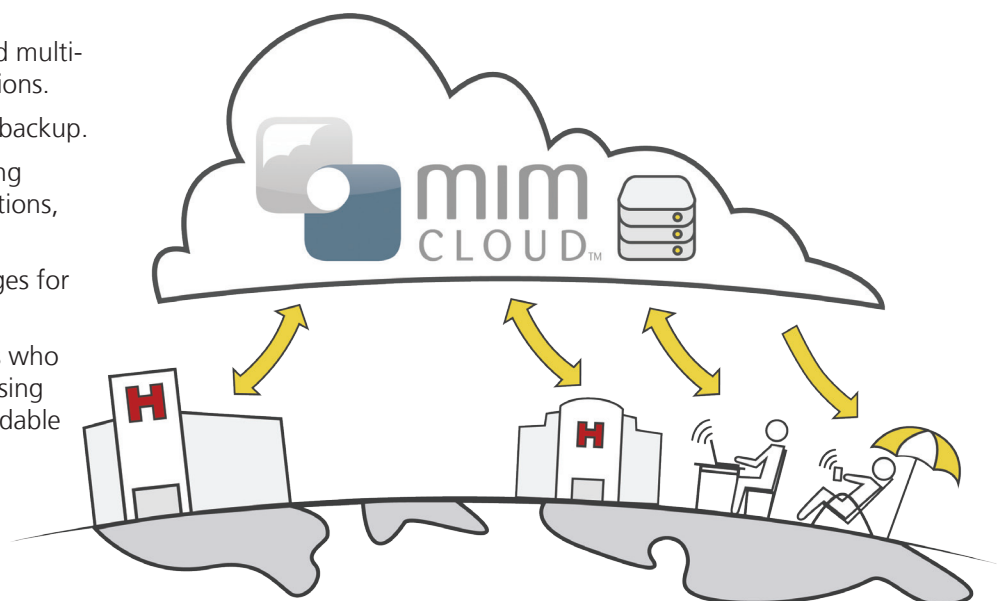

Mobile MIM

- multi-modality fusion & display (CT, MRI, PET, SPECT)
- multi-planar reconstruction of data sets
- multi-touch interface including
 - zoom, fusion, blending, and window/level
- create subtraction/difference images
- data encryption ensures patient privacy
- image access through MIMcloud or a MIM workstation
- available on the iPad, iPhone, and iPod touch

MIMcloud

- access images from anywhere
- web-based viewing solution
- seamless integration with
 - MIM workstation software and Mobile MIM
- share easily with colleagues or patients
- short-term or long-term remote storage
- HIPAA compliant

Apple, the Apple logo, iPad, iPhone, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

MIMpacs / Networking

- DICOM send/receive & query/retrieve
- RT DICOM Support (Plan, Dose, Beam, etc.)
- FTP send/receive
- send volumes and secondary captures to PACS
- fully integrated multi-modality PACS
- secure network user login (LDAP integration)
- automatic study routing / retrieval / archival
- per-user configuration

MIM Workstation / Server System Recommendations - PC

- Intel® Core i7 (Quad Core)
- 8+ GB 800+ MHz RAM
- 16x CD/DVD +/-RW Drive
- 2-8 TB Hard Drive
- 512 MB Dual DVI Graphics Card
- One 24" LCD or Two 19" LCD
- Gigabit Ethernet
- Microsoft Windows 7® Professional 64-bit or
- Microsoft Windows Vista® Business 64-bit

MIM Workstation / Server System Recommendations - Mac

- Mac® Pro – Intel Quad-Core Xeon / iMac® 27-inch – Intel® Core i7
- 8+ GB RAM
- SuperDrive®
- 2-8 TB Hard Drive
- 512 MB ATI or NVIDIA® Graphics Card
- Apple® Display(s)
- Mac OS® X 10.6



MIM Software Inc.

25200 Chagrin Blvd., Suite 200
Cleveland, OH 44122

866-421-2536

www.mimsoftware.com
info@mimsoftware.com

© Copyright 2011 MIM Software Inc.

Beijing MIM Software Corp. Ltd.

北京明维视景医疗软件开发有限公司
地址：北京市海淀区苏州街29号
维亚大厦605室，邮编100080

电话 86-10-82626590

邮箱 info@mimsoftware.com



Emergo Europe

Molenstraat 15
2513 BH The Hague, The Netherlands

Tel: (31) (0) 70 345-8570

Fax: (31) (0) 70 346-7299

service@emergogroup.com



Rev. Date: 4-11-2011